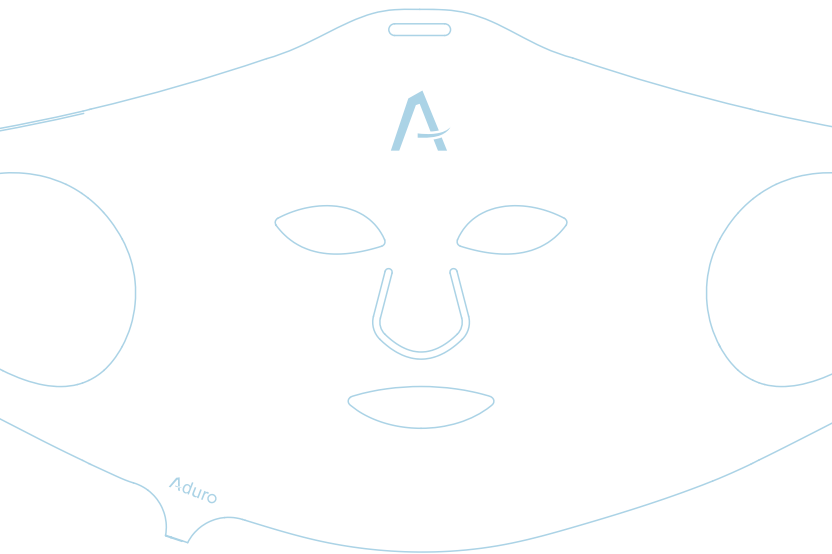
# LED LIGHT THERAPY Facial Mask

User Manual  
Manuel de l'utilisateur  
Bedienungsanleitung  
Manual del usuario  
Manuale utente  
Gebruikershandleiding  
Manual do Utilizador  
Användarhandbok  
用户手册

English	-02
Français	-13
Deutsch	-24
Español	-35
Italiano	-46
Nederlands	-57
Português	-68
Svenska	-79
简体中文	-90

**Aduro**  
English

## 1. The Benefits of Light Therapy

Your skin has the wonderful ability to absorb light energy. Light energy makes skin cells produce more adenosine triphosphate (ATP) stimulating cell renewal, making your skin look younger and brighter. Light therapy is a safe and effective treatment for a variety of skin conditions.

### Aduro 7+1

The Aduro 7+1 LED Facial Mask is designed with 11 different modes, utilizing state-of-the-art infrared technology that boosts vitality and amplifies anti-aging treatments.

#### MODES AND TREATMENTS

WRINKLE REDUCER	ACNE BUSTER	SKIN BALANCER	RED AWAY
<b>Anti-Aging</b> Combats wrinkles and firms the skin.	<b>Anti-Acne</b> Eliminates acne-causing bacteria and reduces acne breakouts.	<b>Evens Skin Tones</b> Helps even skin complexion for naturally beautiful skin.	<b>Reduce Redness</b> Helps alleviate inflammation, sunburns, and rosacea.
RADIANT SKIN	SKIN SOOTHER	CELL REJUVENATOR	SKIN BOOSTER
<b>Revitalization</b> Promotes a brighter glow for more radiant skin and helps relieve rosacea	<b>Soothing</b> Helps to reduce swollen capillaries	<b>Cell Rejuvenation</b> Promotes cell health and regeneration.	<b>Red/Infrared</b> Enhances cellular function, promotes healing, and boosts collagen production.
PAIN BLOCKER	DAILY DOSE		QUICK FIX
<b>BLUE/INFRARED</b> Blue and Infrared light produces NO2 to block pain signals.	<b>RED/BLUE/GREEN/ORANGE/CYAN/PURPLE/YELLOW</b> Combines 7 spectrums of light for full benefits without UV.		<b>RED/ORANGE/GREEN/INFRARED</b> Pulsing of light modes thru Wrinkle Reducer, Revitalization, and Skin Booster for a quick anti-aging lift.

## 2. Company Information

### Company

Aduro Co., LTD.

### Head Office

Spinnerij 67 1185 ZS Amstelveen  
The Netherlands

### Factory

Building A, No. 40, Fuxin Street, Huaide Community  
Bao'an District, Shenzhen, Guangdong Province  
China

### Hong Kong Office

Flat F, 3/F., Hung Wai Industrial Building, 3 Hi Yip Street, Yuen Long, Hong Kong

### USA Address

382 NE 191 Street, Miami, FL 33179 USA

### Support

Email: info@aduroled.com

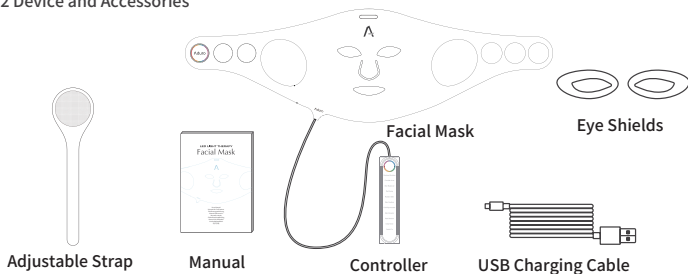
Telephone: +31-6-1473-3907

### 3. Production Information

#### 3.1 Device Specifications

LED color	Blue, Green, Red, Infrared
Output optical power density	30 mw / cm <sup>2</sup>
LEDs	66 Quad LEDs
Operating temperature	0°C ~ 40°C / 32°F ~ 104°F
Storing and transport temperature	-5°C ~ 55°C / 23°F ~ 131°F
Operating humidity	30% ~ 95% without condensation
Storage and Transport Humidity	20% ~ 95% without condensation
Charging supply	5V=2A
Charging port	USB-C
Charging cable	USB-C to USB-A
Charging cable length	101.5cm / 3.33 feet
Battery	2600mAh
Battery type	Lithium-ion
Charging time	240 minutes / 4 Hours
Recommended use	Once a day
Treatment time	10 minutes
Dimensions of the device	706*227mm
Weight of the device	370g

#### 3.2 Device and Accessories

### 3.3 Adverse Reactions

If you have any adverse reactions such as discomfort, skin discoloration, or skin irritation from using this device immediately stop using it and consult a medical professional.

Skin irritation, bruising, contact burns, dizziness, nausea, and headaches are rare but possible side effects. If you experience any of these symptoms immediately stop using the device and contact a medical professional.

### 3.4 Contraindications

- DO NOT use the device if you have cancer.
- DO NOT use the device if you are pregnant or breastfeeding.
- DO NOT let minors operate the device.
- DO NOT use the device on phlebitis, thrombophlebitis, and varicose veins.
- DO NOT use the device if you have a known allergic reaction to silicon. The device contains plastics and silicon that come into contact with your skin during use. If you have any uncertainties consult a medical professional before use.
- DO NOT use the device if you have a photosensitivity disorder (sensitivity to light).
- DO NOT use the device if you are taking any medications that cause photosensitivity. If you have any uncertainties consult a medical professional before use.
- DO NOT use the device if you suffer from light-induced headaches.
- DO NOT use the device if you have a medical history of seizures or epilepsy triggered by light.

Note: For reasons of hygiene, we do not recommend sharing your device with others.

Note: We highly suggest using the device with caution if you are taking any medications. If you have any uncertainties, please consult a medical professional.

Note: If you are currently under the constant care of a physician, consult your physician before using this device.

### 3.5 Warnings

- STOP using the device and consult a medical professional if you experience any adverse reactions from using the device.
- STOP using the device if you experience symptoms of hypotension (dizziness or light-headedness).
- STOP using the device if it is overheating. Turn it off immediately.
- STOP the device if it feels hot or your skin feels uncomfortable.
- DO NOT let children and minors operate the device.

- DO NOT use the device to treat rosacea, moles, warts, open sores, or cancerous lesions.
- DO NOT use the device on areas of your skin without pigmentation or areas that have vitiligo.
- DO NOT shine the light from the LEDs directly into your eyes. This may cause temporary loss of sight.
- DO NOT place the device on any open wounds.
- DO NOT use the device outdoors.
- DO NOT sit, lay down, or step on the device.
- DO NOT drop or insert any objects into any openings on the device.
- DO NOT use the device if it is not working properly or has been dropped or damaged. For warranty and repairs contact customer support.
- DO NOT use or store the device on hot surfaces.
- DO NOT charge the device close to water.
- DO NOT place the device in direct sunlight.
- DO NOT attempt to modify the device in any way. Modifying or opening the device will void your warranty. The device contains no serviceable parts. The batteries are not replaceable.
- DO NOT perform maintenance on or take apart the device. In the case of a malfunctioning device contact customer support.
- DO NOT charge the device if the charging cable is damaged.
- DO NOT use the device with 3rd party accessories.
- DO NOT use solvents or abrasive cleaners on the device.
- DO NOT use the device in any other way than what is specified in this manual.
- DO NOT use the device if you recently had a surgical procedure.
- DO NOT use the device over a bacterial infection as this may result in it worsening.
- DO NOT use the device on painful areas of your body. Consult a medical professional if you experience pain during use.
- DO NOT use the device with a 3rd party charging cable.
- DO NOT use the device in combination with any topicals not endorsed by Aduro.

Note: The device should only be supplied with 5V DC low voltage electricity. If any part of the charging cable is damaged, discontinue use and contact customer support.

### 3.6 Precautions

- Do not scratch or damage the light-emitting surface.
- Not operating the device in accordance with the instructions in the manual will void the warranty.
- Read all instructions thoroughly before use.
- Store the device according to the instructions in the manual.
- This device is intended to be used by one person at a time.
- This device is designed for cosmetic use only. Please contact a medical professional if you have any uncertainties.
- The battery is not removable.
- If the device is to be disposed of, please dispose of it properly.
- Do not use the device if you see any external signs of damage.
- The device should only be powered with 5V DC low voltage electricity. If any part of the charging cable is damaged, discontinue use.
- Consult a medical professional before using the device if you have a fever, acute inflammation, a severe circulatory disorder, light sensitivity, are anesthetized, or have any desensitized areas on your body. If damaged, discontinue use.
- We highly suggest taking caution when using our device while taking any medications. If you are uncertain, consult with a medical professional."

### 3.7 Travel and International Use

The device's lightweight design makes it a perfect travel companion. Here are some travel tips for you:

- If you are travelling to a different region from where you bought the device, please make sure that the device is only supplied with 5V DC voltage electricity.
- Confirm with the flight crew that the device can be carried and/or used on the airplane prior to flight.

### 3.8 Disclaimer

DO NOT use this device to treat any conditions not specified in this manual.

This device is not intended for use in the diagnosis, prevention, or cure of any disease or other medical condition.

## 4. Risks

### 4.1 Potential Risks

Immediately after treatment your skin may appear red. This is normal and the redness will fade within 1-2 hours.

You may experience a temporary mild headache after light treatment. If the headache persists contact a medical professional.

## 5. Light Sensitivity Test

We highly recommend users perform a light sensitivity test before their first use of the device.

1. Place the device on your inner forearm.
2. Turn it on.
3. Wait 3 minutes until the session is finished.
4. If your skin feels hot or uncomfortable during the treatment immediately stop! You might be sensitive to light.
5. The device will automatically turn off after the session is finished. Remove the device and wait for at least 24 hours.
6. Check your skin regularly over 24 hours, if your skin appears red and/or itchy discontinue use of the device.

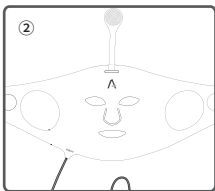
In the case of light sensitivity please contact customer support.

## 6. HOW TO USE THE DEVICE


### 6.1 General Use


Before starting the treatment, thoroughly cleanse your skin, and remove all make up. Not doing this may affect light output and decrease the effectiveness of the treatment.

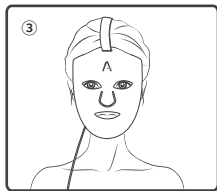Feed the adjustable strap through the hole on the top of the mask ensuring the Velcro can be stuck onto the back of the strap.

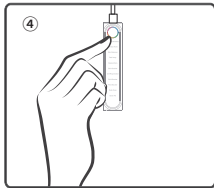You are now ready for your treatment. Put on the mask and adjust the straps to secure it to your face, then adjust the top strap to your desired tightness.

Selectable modes are chosen via the controller. Push the mode button to cycle through the different light modes.

After your light therapy treatment, moisturize and hydrate your skin to lock in the benefits.

## 6.2 Treatment Protocol

The Aduro Facial Mask is designed for a maximum use of 10 minutes daily.

### What to expect during a treatment:

Users may feel a warm sensation on the skin, with increased hydration and plumping from the first treatment.

### After 1 week of treatments:

Most users can expect to notice an improvement in skin texture, with their skin appearing more healthy and smooth.

### After 2 weeks of treatments:

Most users can expect better skin conditions when compared to before treatment started.

### After 4 weeks of treatments:

Many users will notice a more significant skin improvement compared to the weeks prior, with reduced fine lines and wrinkles, acne recovery, and better healthier-looking skin.

## 6.3 Charging

When charging the device, please plug the USB-C end of the data cable into the controller box and the USB-A end into the USB power supply.

When the device is charging, the LED indicator on the controller box will flash red, and when the device is fully charged the indicator will stay green.

Before charging make sure that the plug and wall socket are completely dry. Failure to do so may result in electric shock, short circuit, or fire.

Note: Always unplug the device before use, the device will not work while charging.

## 6.4 Care and Storage

- The lifetime of the lithium battery can be maximized by allowing it to fully charge and then fully discharge during regular use.
- Store the device at room temperature in a dry environment.
- Store the device in its original packaging. Avoid bending, stretching or crushing the device while it is stored or when traveling with it.
- Avoid storing the device in a dusty environment, under direct sunlight, or in damp and humid conditions.
- The device contains no serviceable parts.
- Clean the device with it switched OFF before each use following the cleaning instructions.

## 6.5 Cleaning

- Clean the device with a water-based solution or damp cleaning cloth. Stubborn stains can be removed by wiping gently with 70% alcohol.
- DO NOT use any cleaning appliances on the device as these may damage it.
- DO NOT clean the device while it is turned ON.
- DO NOT clean the device while charging.

## 6.6 Troubleshooting

If the solutions suggested in this section fail to fix the problem, and your device is within the warranty period, please contact customer support.

△ DO NOT attempt to diagnose the device yourself.

△ DO NOT attempt to modify the device in any way. This device has NO serviceable parts.

Problem	Solution
The device does not switch ON.	The device is not fully charged, because the charging cable was not properly inserted into the device. Insert the charging cable properly into the device and then charge it.
	The device is not fully charged, because the charging cable was not properly inserted into it. Insert the charging cable properly into the device and then charge it.
	The power outlet does not have power. Try a different power outlet.
	The ON/OFF button has not been pushed correctly. Hold down the ON/OFF button.
The device switches OFF unexpectedly.	The battery level is too low. Charge the battery.

## 6.7 Disposal

DO NOT dispose of the device together with household waste. The device should be properly disposed of because it contains a lithium-ion battery. For protection of the environment, please safely dispose of the device following your local regulations.

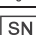
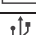
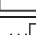
## 7. Warranty

Aduro guarantees that this product will be free from material and quality defects within one year from the date of purchase. This product is not warranted as follows:

- This product's warranty does not cover damage caused by misuse or abuse.
- This product's warranty does not cover damage caused by accident.
- This product's warranty does not cover damage caused by the attachment of any unauthorized or third-party accessories.
- This product's warranty does not cover damage caused by alteration to the product.
- This product's warranty does not cover damage caused by the user after purchase.
- This warranty does not cover broken light-emitting diodes unless the number exceeds 4.
- Aduro will not be responsible for any incidental, consequential, or circumstantial damage.
- To obtain warranty service for this product, please return the product (postage paid by you) to Aduro with a sales receipt as proof of purchase. Once the item has been received, Aduro will repair or replace it following our warranty policy, then return your product to you (postage included). Warranty is only valid through Aduro approved service centers.
- This warranty does not affect your statutory rights.

## 8. Symbols

Symbol	Explanation
	CE certified and authorized to be sold in Europe.
	Electrical and Electronic Equipment Waste (WEEE), indicates the device may not be discarded as household waste (2012/19/EU).
	Electromagnetic radiation from the device is below the limits specified by the Federal Communications Commission, and the manufacturer has followed the requirements of the Supplier's Declaration of Conformity authorization procedures.

	Consult instructions for use.
	Conformity to Australian Directives on health, safety and environmental protection standards. RCM stands for Regulatory Compliance Mark.
	This symbol means “Manufactured by” and includes the name and address of the legal manufacturer.
	Date of manufacture.
	Keep the device dry.
	Serial number.
	Catalogue number.
	Batch code.
	USB connector.
	The device uses a lithium-ion polymer battery.
	Made in China.
	Keep the device away from direct sunlight.
	Safe operating temperatures.
	Safe operating humidity.
	The cardboard packaging that the device came in can be recycled.
	Type BF Applied Part, B indicates it is used for the body and F indicates that it is not connected to the mains voltage.
	Unique Device Identifier.

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Spinnerij 67 1185 ZS Amstelveen  
The Netherlands

### Usine

Building A, No. 40, Fuxin Street, Huaide Community  
Bao' an District, Shenzhen, Guangdong Province  
China

### Bureau de Hong Kong

Flat F, 3/F., Hung Wai Industrial Building, 3 Hi Yip Street, Yuen Long, Hong Kong

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