
LED LIGHT THERAPY

Palm

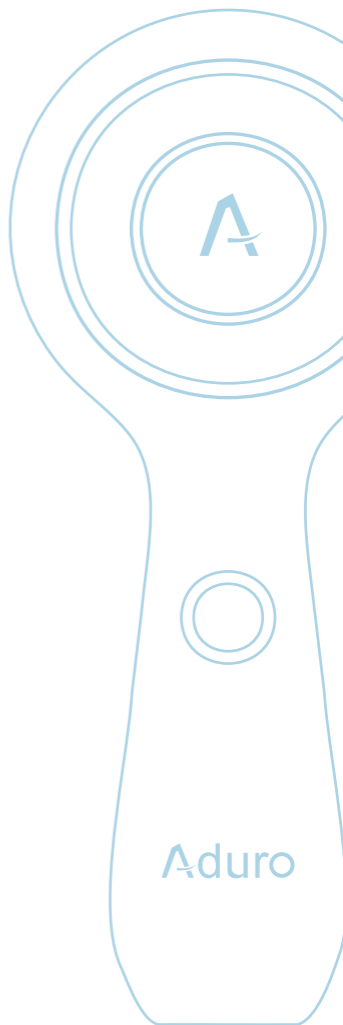User Manual

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Manual do Utilizador

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English	-02
Français	-10
Deutsch	-18
Español	-26
Italiano	-34
Nederlands	-42
Português	-50
Svenska	-58
简体中文	-66

Aduro
English

1. Company Information

Company

Aduro Co., LTD.

Head Office

Spinnerij 67 1185 ZS Amstelveen
The Netherlands

Factory

Building A, No. 40, Fuxin Street, Huaide Community
Bao'an District, Shenzhen, Guangdong Province
China

Hong Kong Office

Flat F, 3/F., Hung Wai Industrial Building, 3 Hi Yip Street, Yuen Long, Hong Kong

USA Address

382 NE 191 Street, Miami, FL 33179 USA

Support

Email: info@aduroled.com

Telephone: +31-6-1473-3907

2. Product Information

2.1 What is the Aduro PALM used for?

The Aduro PALM is designed to make skin imperfections disappear, pioneering a new way of using light therapy with the combination of cosmetics to improve the appearance and vitality of your skin.

Aduro Mode 1: Red Light+ Infrared Light

This is designed to combat signs of aging, such as wrinkles and fine lines, through the use of red&infrared light technology.

Aduro Mode 2: Red Light+ Blue Light

This is designed to combat acne, helping speed up recovery and the healing time of acne breakouts, through the use of red&blue light technology.


2.2 Device Specifications

LED color	Red, Infrared, Blue
Output optical power density	60 mw/cm ²
LEDs	12 Triple LEDs
Operating temperature	0°C to 40°C 32°F to 104°F
Storing and transport temperature	-5°C to 55°C 23°F to 131°F
Operating humidity	30% to 95% without condensation
Storage and Transport Humidity	20% to 95% without condensation
Charging supply	5V =2A
Charging port	USB-C
Charging cable	USB-C to USB-A
Charging cable length	1.2 meter / 3.937 feet
Battery	3.7V, 2600mAh
Battery type	Lithium-ion
Charging time	6.5h
Recommended use	Once a day
Treatment time	3 minutes
Dimensions of the device	173.28 x 74.4 mm
Weight of the device	133g

2.3 Device Accessories

Palm

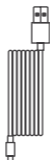USB Charging Cable

Manual

2.4 Adverse Reactions

- If you have any adverse reactions such as discomfort, skin discoloration, or skin irritation from using this device immediately stop using it and consult a medical professional.
- Skin irritation, bruising, contact burns, dizziness, nausea, and headaches are rare but possible side effects. If you experience any of these symptoms immediately stop using the device and contact a medical professional.

2.5 Contraindications

- DO NOT use the device if you have cancer.
- DO NOT use the device if you are pregnant or breastfeeding.
- DO NOT let minors operate the device.
- DO NOT use the device on phlebitis, thrombophlebitis, and varicose veins.
- DO NOT use the device if you have a known allergic reaction to silicon. The device contains plastics and silicon that come into contact with your skin during use. If you have any uncertainties consult a medical professional before use.
- DO NOT use the device if you have a photosensitivity disorder (sensitivity to light).
- DO NOT use the device if you are taking any medications that cause photosensitivity. If you have any uncertainties consult a medical professional before use.
- DO NOT use the device if you suffer from light-induced headaches.
- DO NOT use the device if you have a medical history of seizures or epilepsy triggered by light.

Note: For reasons of hygiene, we do not recommend sharing your device with others.

Note: We highly suggest using the device with caution if you are taking any medications. If you have any uncertainties, please consult a medical professional.

Note: If you are currently under the constant care of a physician, consult your physician before using this device.

2.6 Warnings

- STOP using the device and consult a medical professional if you experience any adverse reactions from using the device.
- STOP using the device if you experience symptoms of hypotension (dizziness or light-headedness).
- STOP using the device if it is overheating. Turn it off immediately.
- STOP the device if it feels hot or your skin feels uncomfortable.

- DO NOT let children and minors operate the device.
- DO NOT use the device to treat rosacea, moles, warts, open sores, or cancerous lesions.
- DO NOT use the device on areas of your skin without pigmentation or areas that have vitiligo.
- DO NOT shine the light from the LEDs directly into your eyes. This may cause temporary loss of sight.
- DO NOT place the device on any open wounds.
- DO NOT use the device outdoors.
- DO NOT sit, lay down, or step on the device.
- DO NOT drop or insert any objects into any openings on the device.
- DO NOT use the device if it is not working properly or has been dropped or damaged. For warranty and repairs contact customer support.
- DO NOT use or store the device on hot surfaces.
- DO NOT charge the device close to water.
- DO NOT submerge the device in water.
- DO NOT place the device in direct sunlight.
- DO NOT attempt to modify the device in any way. Modifying or opening the device will void your warranty. The device contains no serviceable parts. The batteries are not replaceable.
- DO NOT perform maintenance on or take apart the device. In the case of a malfunctioning device contact customer support.
- DO NOT charge the device if the charging cable is damaged.
- DO NOT use the device with 3rd party accessories.
- DO NOT use solvents or abrasive cleaners on the device.
- DO NOT use the device in any other way than what is specified in this manual.
- DO NOT use the device if you recently had a surgical procedure.
- DO NOT use the device over a bacterial infection as this may result in it worsening.
- DO NOT use the device on painful areas of your body. Consult a medical professional if you experience pain during use.
- DO NOT use the device in combination with any topicals not endorsed by Aduro.

Note: The device should only be supplied with 5V DC low voltage electricity. If any part of the charging cable is damaged, discontinue use and contact customer support.

2.7 Precautions

- Do not scratch or damage the light-emitting surface.
- Not operating the device in accordance with the instructions in the manual will void the warranty.
- Read all instructions thoroughly before use.
- Store the device according to the instructions in the manual.
- This device is intended to be used by one person at a time.
- This device is designed for cosmetic use only. Please contact a medical professional if you have any uncertainties.
- The battery is not removable.
- If the device is to be disposed of, please dispose of it properly.
- Do not use the device if you see any external signs of damage.
- The device should only be powered with 5V DC low voltage electricity. If any part of the charging cable is damaged, discontinue use.
- Consult a medical professional before using the device if you have a fever, acute inflammation, a severe circulatory disorder, light sensitivity, are anesthetized, or have any desensitized areas on your body.
- We highly suggest taking caution when using our device while taking any medications. If you are uncertain, consult with a medical professional.

2.8 Travel and International Use

The device's lightweight design makes it a perfect travel companion. Here are some travel tips for you:

- If you are travelling to a different region from where you bought the device, please make sure that the device is only supplied with 5V DC low voltage electricity.
- Confirm with the flight crew that the device can be carried and/or used on the airplane prior to flight.

2.9 Disclaimer

DO NOT use this device to treat any conditions not specified in this manual.

This device is not intended for use in the diagnosis, prevention, or cure of any disease or other medical condition.

3. Risks

3.1 Potential Risks

Immediately after treatment your skin may appear red. This is normal and the redness will fade within 1-2 hours.

You may experience a temporary mild headache after light treatment. If the headache persists contact a medical professional.

4. Light Sensitivity Test

We highly recommend users perform a light sensitivity test before their first use of the device.

1. Place the device on your inner forearm.
2. Turn it on.
3. Wait 3 minutes until the session is finished.
4. If your skin feels hot or uncomfortable during the treatment immediately stop! You might be sensitive to light.
5. The device will automatically turn off after the session is finished. Remove the device and wait for at least 24 hours.
6. Check your skin regularly over 24 hours, if your skin appears red and/or itchy discontinue use of the device. In the case of light sensitivity please contact customer support.


5. How To Use The Device

5.1 General Use

① Thoroughly clean the area of the skin that you want to treat, hold down the power button for 2 seconds to start. Quickly press the power button to switch treatment modes.

② Place the PALM over the area to be treated and hold close to the skin. One session will last 3 minutes.

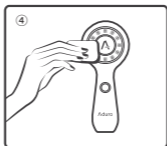③ After finishing one area, you can move to the next area to be treated and repeat on as many treatment areas as you wish.

④ Clean the device with a water-based solution or damp cleaning cloth. Ensure the device is NOT connected to the power supply when cleaning.

5.2 Treatment Protocol

Use the device 1 time per day for a minimum of 4 weeks. After that you can continue using the device as required based upon changes in your skin condition. It is not recommended to use the device more than 7 times per week.

5.3 Charging

Avoid fast charging as this can damage the battery.

Before charging make sure that the plug and wall socket are completely dry. Failure to do so may result in electric shock, short circuit, or fire.

Note: Always unplug the device before use, the device will not work while charging.

5.4 Care and Storage

- The lifetime of the lithium battery can be maximized by allowing it to fully charge and then fully discharge during regular use.
- Store the device at room temperature in a dry environment.
- Store the device in its original packaging. Avoid bending, stretching or crushing the device while it is stored or when traveling with it.
- Avoid storing the device in a dusty environment, under direct sunlight, or in damp and humid conditions.
- The device contains no serviceable parts.
- Clean the device with it switched OFF before each use following the cleaning instructions.

IMPORTANT

The product has an average lifecycle of two years under normal operation. Following the instructions in this manual is considered normal operation.

Handle the device with care.

5.5 Cleaning

- Clean the device with a water-based solution or damp cleaning cloth. Stubborn stains can be removed by wiping gently with 70% alcohol.
- DO NOT use any cleaning appliances on the device as these may damage it.
- DO NOT clean the device while it is turned ON.
- DO NOT clean the device while charging.

5.6 Troubleshooting

If the solutions suggested in this section fail to fix the problem, and your device is within the warranty period, please contact customer support.

△ DO NOT attempt to diagnose the device yourself.

△ DO NOT attempt to modify the device in any way. This device has NO serviceable parts.

Problem	Solution
The device does not switch ON.	The device is not fully charged, because the charging cable was not properly inserted into the device. Insert the charging cable properly into the device and then charge it.
	The device is not fully charged, because the charging cable was not properly inserted into it. Insert the charging cable properly into the device and then charge it.
	The ON/OFF button has not been pushed correctly. Hold down the ON/OFF button.
The device switches OFF unexpectedly.	The battery level is too low. Charge the battery.

5.7 Disposal

DO NOT dispose of the device together with household waste. The device should be properly disposed of because it contains a lithium-ion battery. For protection of the environment, please safely dispose of the device following your local regulations.

6. Warranty

Our devices are warranted to be free from defects in material and workmanship for 12 months from its original purchase date under normal operation. This warranty extends only to the original retail purchaser with original proof of purchase and only when purchased from an authorized retailer or reseller. The warranty coverage is extended if the device is used in conjunction with authorized device accessories. If the device should become defective within the warranty period, contact the retailer that you purchased the device from. This warranty does not cover devices damaged by:

- Accident, misuse, abuse, modification.
- Use with unauthorized 3rd party accessories or other than as indicated.
- Connecting to incorrect current and voltage.

Aduro shall not be responsible for any incidental, special or consequential damages resulting from the use of its devices. All implied warranties (including but not limited to implied warranties of fitness and merchantability) are limited in duration to one year from date of original purchase.

This warranty gives you specific legal rights, and you may also have other rights that vary from country to country. Some countries do not allow the exclusion or limitation of incidental or consequential damages, or limitations on how long an implied warranty lasts, so some of the above limitations and exclusions may not apply to your country.

7. How To Use The Device

Symbol	Explanation
	CE certified and authorized to be sold in Europe.
	Electrical and Electronic Equipment Waste (WEEE), indicates the device may not be discarded as household waste (2012/19/EU).
	Read the enclosed documents before use and take all safety precautions required.
	Consult instructions for use.
	Conformity to Australian Directives on health, safety and environmental protection standards. RCM stands for Regulatory Compliance Mark.
	This symbol means “Manufactured by” and includes the name and address of the legal manufacturer.
	Date of manufacture.
	Keep the device dry.
	Serial number.

	Catalogue number.
	Batch code.
	USB connector.
	The device uses a lithium-ion polymer battery.
	Made in China.
	Keep the device away from direct sunlight.
	Safe operating temperatures.
	Safe operating humidity.
	The cardboard packaging that the device came in can be recycled.
	Type BF Applied Part, B indicates it is used for the body and F indicates that it is not connected to the mains voltage.
	Unique Device Identifier.

MADE IN CHINA

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Entreprise

Aduro Co., LTD.

Siège social

Spinnerij 67 1185 ZS Amstelveen
The Netherlands

Usine

Building A, No. 40, Fuxin Street, Huaide Community
Bao'an District, Shenzhen, Guangdong Province
China

Bureau de Hong Kong

Flat F, 3/F., Hung Wai Industrial Building, 3 Hi Yip Street, Yuen Long, Hong Kong

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